
GET TO KNOW YOUR POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT

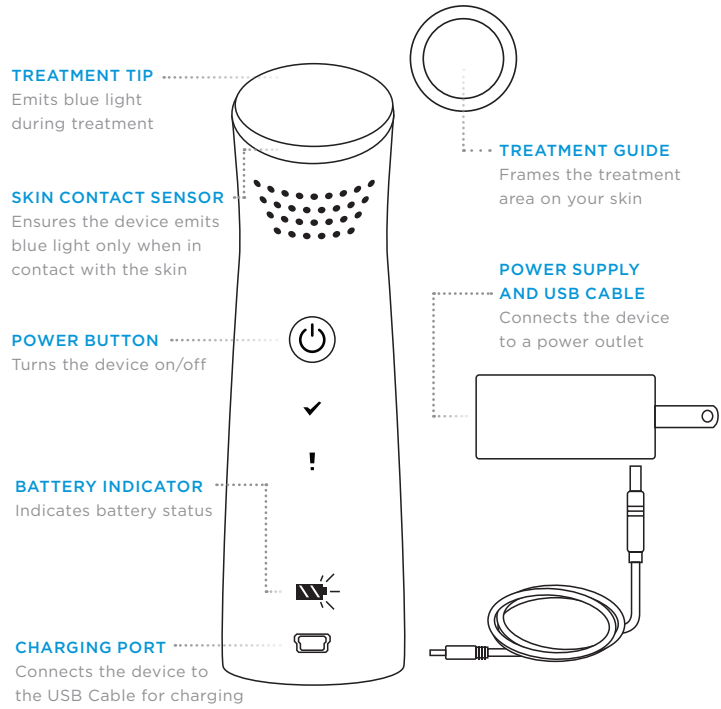tria[®]
**POSITIVELY CLEAR
 ACNE CLEARING BLUE LIGHT**

Quick Start Guide

IMPORTANT DEVICE INDICATORS

TREATMENT COMPLETION INDICATOR
 Lights when treatment of the selected area is complete (5 mins). The device will also chime and then automatically turn off.

DEVICE ERROR INDICATOR
 Lights if there is a device error. The device will also make a buzzing tone.

LOW BATTERY INDICATOR
 Battery indicator flashes. The battery needs to be recharged. When the device is plugged in and charging, the battery indicator will show the battery filling up.

FULL BATTERY INDICATOR
 Battery indicator shows a full battery (3 bars). The device is fully charged and ready to use.

Tria Beauty, Inc.
 4160 Dublin Blvd., Suite 200
 Dublin, CA 94568 USA

TRIA BEAUTY CUSTOMER CARE
 877-321-8742
 customercare@triabeauty.com
 triabeauty.com

tria.

Please refer to the Instructions for Use booklet for complete information about the Positively Clear Acne Clearing Blue Light, including warnings and indications for use.

©2016 Tria Beauty, Inc. All rights reserved. Tria, Tria Beauty and the Tria logo are registered trademarks of Tria Beauty, Inc. Patents pending.

15611B

GET STARTED

BEFORE FIRST TREATMENT

For detailed directions and important safety warnings, please read the complete *Instructions for Use*.

CHARGE


Charge the Positively Clear Acne Clearing Blue Light until the battery indicator stops showing the battery filling up and displays a full battery (less than 6 hrs). Your device may arrive fully charged and not require additional charging before first use.

HOW TO TREAT

Treat once a day and for best results, do not skip treatments. Start by cleansing and drying the skin you wish to treat.

1 APPLY THE TREATMENT GUIDE

Take one treatment guide, and place it on the area of skin that you wish to treat so that it "frames" an area with pimples or pustules. The area in the middle of the guide is the treatment area.

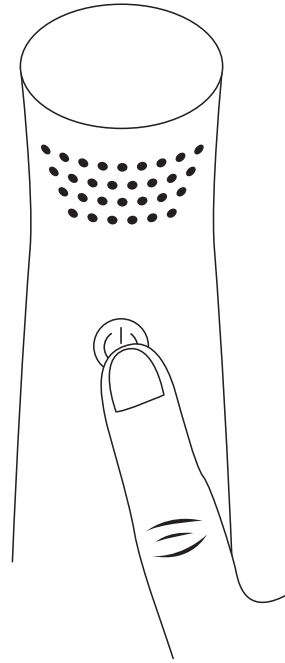2 TURN ON YOUR DEVICE

Before placing the device on the area you wish to treat, press and release the power button. The device will beep and is now ready to use.

NOTE: The device will not turn on while plugged in and will not treat if turned on after it has already been placed on your skin. Make sure to unplug the device and turn on prior to making contact with your skin.

3 START YOUR TREATMENT


Place the **treatment tip against your skin** inside the treatment guide.

Slide the treatment tip continuously and evenly across the entire treatment area, maintaining constant contact with the skin, for 5 minutes.


You will hear a fan whenever the device is emitting blue light treatment.

When treatment of the selected area is complete, the device will chime, the treatment completion indicator (✓) will light and the device will then automatically turn off.

Repeat on other areas as needed.

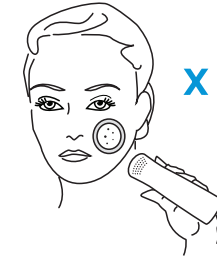X Incorrect skin contact

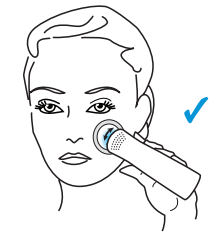✓ Correct skin contact & motion

NOTE: If you lift the device off your skin, the device will stop emitting blue light and the fan will turn off. Place the treatment tip back onto the treatment area to continue treating.

GIVE YOUR ACNE-FIGHTING RESULTS A LIFT

Enhance your Positively Clear Acne Clearing Blue Light experience with our Positively Clear skincare, expertly formulated for use with blue light treatments.

For more information on Tria's acne-fighting skincare products, visit tribeauty.com.

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