tria®

**AGE-DEFYING
EYE WRINKLE CORRECTING LASER**

Quick Start Guide

GET TO KNOW YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER

Age-Defying Eye Wrinkle Correcting Laser

IMPORTANT DEVICE INDICATORS

WHEN THE INDICATOR LIGHT IS . . .	IT MEANS . . .
Green	Steady green: Device is ready to use.
	Slowly flashing green: Device needs charging.
	Rapidly flashing green: Device is charging.
Purple (with vibration)	Device is gliding on skin and emitting laser treatment light.
Blue	Device is locked and must be unlocked for use.
No Light	Treatment is complete. Device has turned off.
Red	Device has an internal error. See troubleshooting guide.

NOTE: You must unplug the device to perform a treatment.

Tria Beauty, Inc.
4160 Dublin Blvd., Suite 200
Dublin, CA 94568 USA

CUSTOMER CARE
877-321-8742
customer care@triabeauty.com
triabeauty.com

tria.

Please refer to the *Instructions for Use* booklet for complete information about the Age-Defying Eye Wrinkle Correcting Laser, including warnings and indications for use.

©2015 Tria Beauty, Inc. All rights reserved. Tria, Tria Beauty and the Tria logo are registered trademarks of Tria Beauty, Inc. Patents pending

14509D

BEFORE FIRST TREATMENT

For detailed directions and important safety warnings, please read the complete *Instructions for Use*.

HOW TO TREAT

Start by cleansing and drying the eye contour area. After treatment, apply a nourishing eye cream.

For best results, use the Tria Priming Cleanser and Nourishing Eye Renewal Cream, which are expertly formulated to complement your nightly laser treatments.


1 TURN ON DEVICE

Press and release the power button. The indicator light will turn green to indicate that the device is ready to use.

2 START YOUR TREATMENT

- Identify the treatment area around your eye as shown.

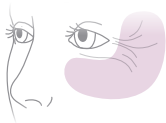
CHARGE: Charge the Age-Defying Eye Wrinkle Correcting Laser until the battery light turns steady green, indicating the device is fully charged (up to 2 hrs). Your device may arrive fully charged and not require additional charging before first use.

START YOUR TREATMENT (CONT.)

- Cover the entire area by placing the laser tip flat against the skin and slowly gliding in small, circular motions. The indicator light glows purple and the device vibrates to indicate you are maintaining good skin contact and treating correctly.
- Treatment will complete after about 1 minute. The indicator light will turn off, the device will vibrate on and off several times, and automatically turn off.
- Repeat on your other eye area.
- Treat 7 days a week for 8 weeks.

 **DO NOT** treat each eye area more than once per day (once every 24 hours). Doing this could result in increased risk of skin irritation.

Circular, Gliding Motion

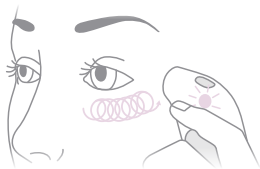 **DO NOT** use to treat your eye or eyelids or aim the device directly into your eyes. Doing this could cause serious eye injury such as a permanent or temporary decrease in vision or permanent blindness. Treat around the eyes only where you can feel bone beneath the skin.

HELPFUL TIPS FOR BEST RESULTS

CLEAN BEFORE USE

Your device is most effective when clean. Ensure optimal results by wiping the laser treatment tip with a dry cotton swab or microfiber cloth prior to each use.

SOOTHE AFTER TREATMENT AND KEEP SKIN HYDRATED

Apply a post-treatment eye cream to help soothe and nourish skin after treatment and apply a moisturizer as needed to keep skin well hydrated.

PRACTICE IN A MIRROR

Perfect your technique by using your device in front of a mirror. The indicator light will glow purple when you are making good skin contact and treating correctly.

TREAT AT NIGHT

Treat in the evening before bedtime to allow skin renewal to begin and side effects to diminish overnight.

PROTECT YOUR SKIN

Protect your skin by applying a broad spectrum sunscreen of at least SPF 30 every morning. Daily sun protection will protect healthy skin created by your Age-Defying Eye Wrinkle Correcting Laser to maintain results longer and help prevent new photodamage.

GIVE YOUR ANTI-AGING RESULTS A LIFT

Enhance your Age-Defying Eye Wrinkle Correcting Laser experience with our cleanser and eye cream, expertly formulated for laser-treated skin.

PRIMING CLEANSER

Prep your skin for treatment. This gentle cleanser, formulated with Willow Bark to promote exfoliation and Birch Tree Sugars to help hydrate skin, ensures skin is perfectly prepped for laser treatment by removing any makeup, sunscreen and other impurities that can block the laser light.

NOURISHING EYE RENEWAL CREAM

Follow your laser treatment with our eye cream to help soothe and hydrate laser-treated skin while reducing puffiness. Specially formulated with Hyaluronic Acid and Ginger Extract to help hydrate, reduce redness and soothe skin. The metal applicator helps calm and cool skin after laser treatment.

For more information on Tria's anti-aging skincare products, visit [triabeauty.com](https://www.triabeauty.com).

A close-up, soft-focus portrait of a woman with blonde hair pulled back, looking slightly to the left. The background is bright and out of focus. The Tria logo is centered over her face.

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